

Catch a new wave to easy cleaning.

In most spinal punches and rongeurs, tissue and debris can collect and remain trapped between the main body and slider. The Clean Wave Punch and Rongeurs have a slider with a wavelike shaped design.

- Wavelike shape allows cleaning utensils to easily reach through the recesses and openings between the slider and main body for cleaning prior to sterilization.
- Wavelike shape also helps to reduce friction between the slider and main body, providing a smooth function.
- Forward angle on punches provides precise cutting action.
- Thin foot plate on 1mm-6mm punches.
- Straight, up and down bite styles on rongeurs.
- Bone ejector is incorporated into the design of 2mm-6mm punches.
- Made from German stainless steel.
- Quality verified and guaranteed by gSource.


Please inquire about the availability of any size and style not shown on this product sheet.



gS 67.0223 Clean Wave Caspar IVD Rongeur pictured on left.
See back of this page for available Clean Wave Rongeurs.

- gS 69.3081** 8" forward 1mm
- gS 69.3082** 8" forward 2mm, with ejector
- gS 69.3083** 8" forward 3mm, with ejector
- gS 69.3084** 8" forward 4mm, with ejector
- gS 69.3085** 8" forward 5mm, with ejector

Clean Wave Kerrison Punch with thin foot plate Ferris-Smith-Kerrison handle 5"

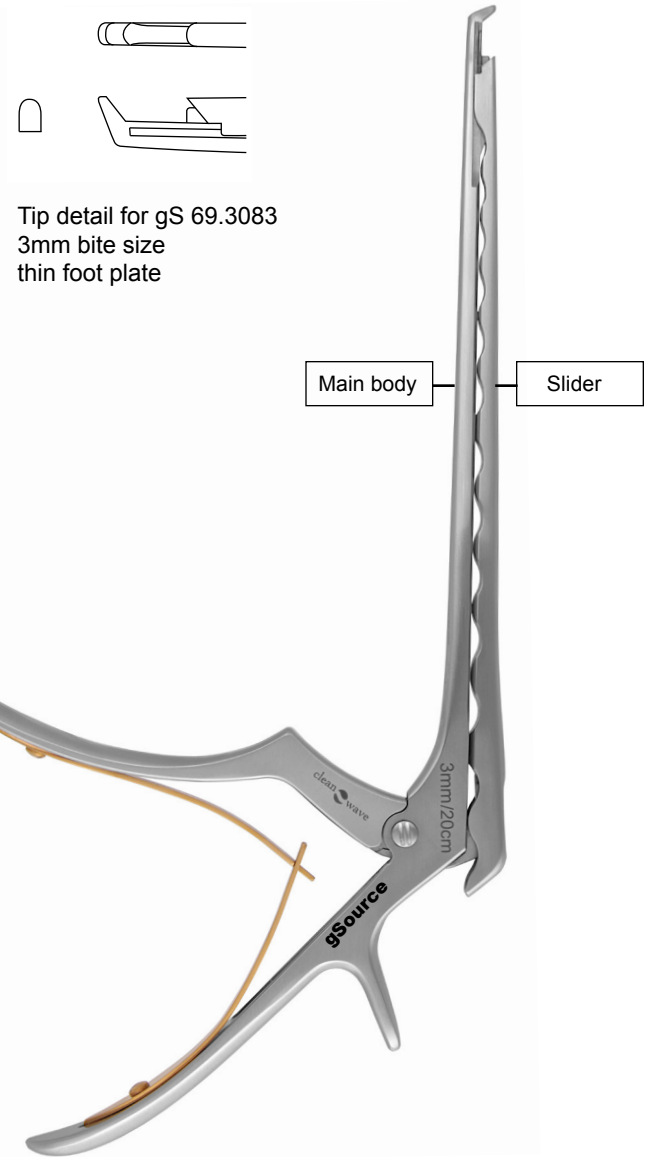

gS 69.3083 Clean Wave Kerrison Punch is pictured above.

See back of this page for available Clean Wave Rongeurs.


Distributed by:

Medical Resources Group Limited
1925 South Main Street
Morton, IL 61550
O: 888-674-6621
F: 309-263-6903
E: sales@mrgona.com



ISO 13485: 2003 Certified
FS 589741
FDA Registered

Jaw bite styles

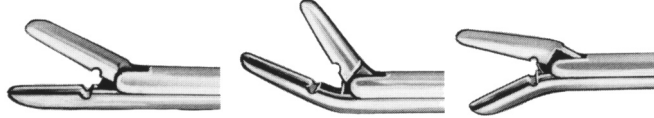

straight

up

down

- gS 67.0222** straight, 2x12mm
- gS 67.0223** straight, 3x12mm
- gS 67.0224** straight, 4x14mm

- gS 67.0242** up, 2x12mm
- gS 67.0243** up, 3x12mm
- gS 67.0244** up, 4x14mm

- gS 67.0232** down, 2x12mm
- gS 67.0233** down, 3x12mm
- gS 67.0234** down, 4x14mm

Clean Wave Caspar IVD Rongeur 7 1/4" fenestrated, serrated side jaws

- gS 67.8300** straight
- gS 67.8320** up
- gS 67.8340** down

Clean Wave Cushing Rongeur 7" 2mm

- gS 67.8800** straight
- gS 67.8820** up
- gS 67.8840** down

Clean Wave Love Gruenwald Rongeur 7" 3mm


- gS 67.9400** straight
- gS 67.9420** up
- gS 67.9440** down

Clean Wave Spurling Rongeur 7" 4mm

Medical Resources Group Limited
1925 South Main
Morton, IL 61550

O: 888-674-6621
F: 309-263-6903
E sales@mrgona.com

www.mrgona.com

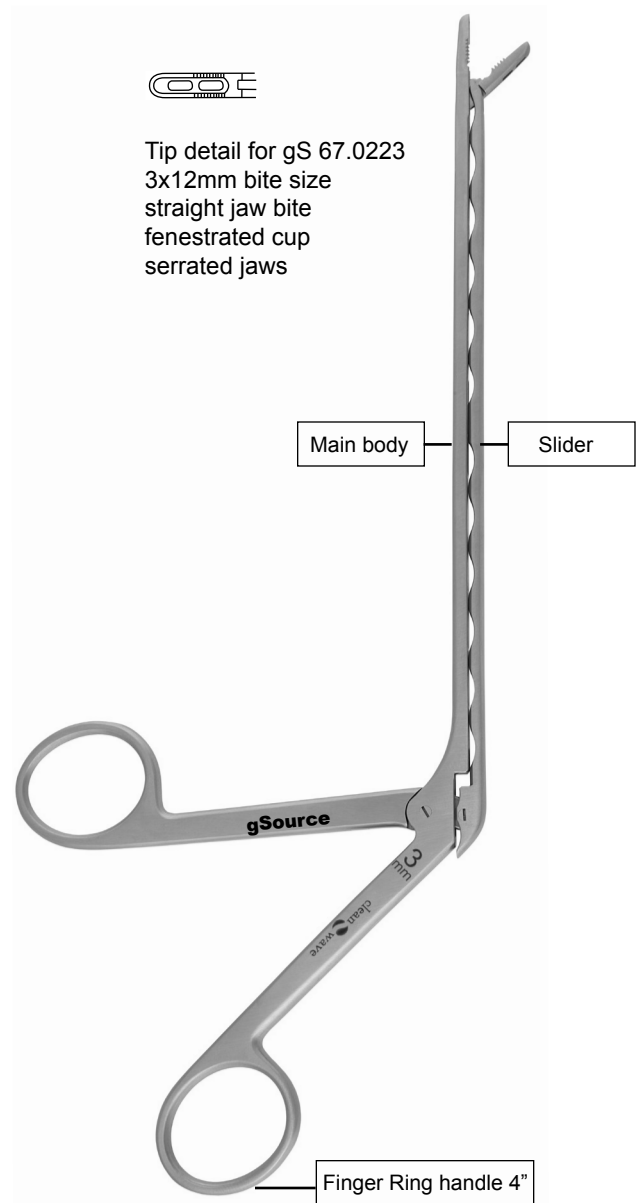

Tip detail for gS 67.0223
3x12mm bite size
straight jaw bite
fenestrated cup
serrated jaws

gS 67.0223 Clean Wave Caspar IVD Rongeur is pictured above.
See front of this page for Clean Wave Kerrison Punches available.



ISO 13485: 2003 Certified
FS 589741
FDA Registered